


Clonality and Lymphomas

Within the oncology field Experteam offers kits to evaluate the polyclonal or monoclonal origin of T or B cell lymphomas. The technique for B cell analysis is based on specific amplification of V-D-J regions within the immunoglobulin heavy chain genes. Similar to the technique used in B cell lymphoma analysis is the technique for T cell lymphoma analysis using amplification of the V-J regions of the gamma chain genes of the T cell receptor.

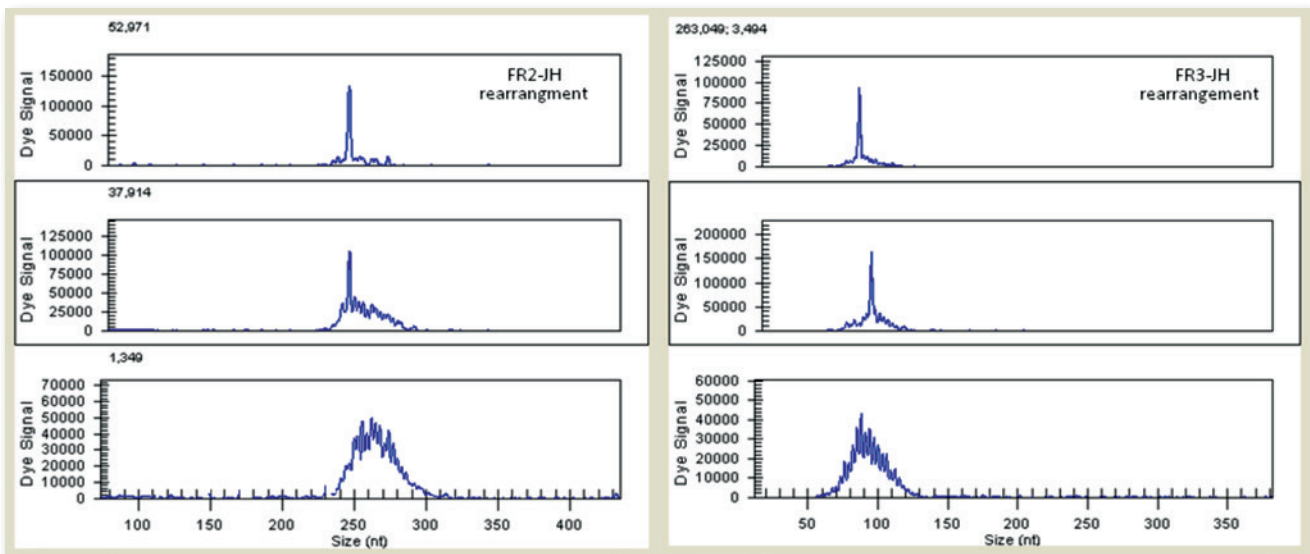

Examples of positive and negative samples for detection of clonal immunoglobulin heavy chain (IgH) gene rearrangements.

Code	Name of kit	Technical specs	Amount of tests	Info	Additional info
BL.01	B cell lymphoma kit	Agarose Gel Electrophoresis Kit	40	CE/IVD	
BL.01FL	B cell lymphoma kit-FL	Capillary Electrophoresis Kit	40	CE/IVD	Appl. Bios. & Beckman Coult.
TL.01	T cell lymphoma kit	Agarose Gel Electrophoresis Kit	40	CE/IVD	
TL.01FL	T cell lymphoma kit-FL	Capillary Electrophoresis Kit	40	CE/IVD	Appl. Bios. & Beckman Coult.

